

# TECHNICAL DATA SHEET



## ANATOMICALLY SHAPED PADS

### SAN SENI



## I. PRODUCT DESCRIPTION

Anatomically shaped pads for people with moderate and heavy incontinence needs. Two anatomically shaped absorbent cores are made of cellulose fluff pulp (in Prima size – one absorbent core) with superabsorbent. Contains distributing layer (EDS). Side gathers with elastic yarn along absorbent core are directed outside preventing side leakages. Wetness indicator – black printed inscription that fades and two yellow hot-melt adhesive stripes change colour into green, when exposed to wetness. Vapour-permeable outer layer on the whole surface.

## II. RAW MATERIALS

Cellulose, superabsorbent, vapour-permeable white laminate with two stripes in the respective colours, hydrophilic-hydrophobic nonwoven, hydrophobic nonwoven, distribution layer (EDS), elastic yarn, ink of indication, hot-melt adhesives, natural green tea extract.

## III. DIMENSIONS\*

Size	Total length [mm]	Length of absorbent core [mm]		Total width [mm]	Width of absorbent core [mm] (front-crotch-back)
		top	bottom		
Prima	610	-	480	325	250-140-260
Normal	625	370	580	320	240-120-260
Uni	625	370	580	320	240-120-260
Regular	690	400	610	370	290-170-330
Maxi	690	400	610	370	290-170-330
Plus	690	400	610	370	290-170-330
Plus Extra	750	450	660	390	310-170-350

\* tolerance of dimensions +/- 10 mm

## IV. PROPERTIES

Size	Weight [g] no less than	Absorbency [g] no less than (ISO 11948-1)	Retention [g] average (NAFC)
Prima	48,8	1100	240
Normal	55,1	1300	410
Uni	81,7	1700	530
Regular	66,5	1500	420
Maxi	110,2	2200	580
Plus	126,3	2800	850
Plus Extra	130,2	3500	1000

## V. BASIC LOGISTIC DATA

Size	Code	Package [pcs]	Packages per box	Boxes per pallet	Box type
Prima	SE-093-PR30-001	30	4	28	carton
Normal	SE-093-NO10-001	10	8	32	carton
	SE-093-NO30-001	30	4	28	carton
Uni	SE-093-UN10-004	10	8	32	carton
	SE-093-UN30-001	30	4	24	carton
Regular	SE-093-RE30-001	30	3	24	carton
Maxi	SE-093-MA10-001	10	8	24	carton
	SE-093-MA30-001	30	3	24	carton
Plus	SE-093-PL10-001	10	8	20	carton
	SE-093-PL30-001	30	3	20	carton
Plus Extra	SE-093-PX30-001	30	3	16	carton

Cartons with pads must be placed on wooden pallet 1200x800 mm, covered by PE sheet (from the bottom and from the top) wrapped together with PE stretch film.

## VI. RECOMMENDED STORAGE CONDITIONS

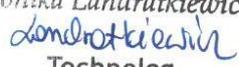
It is recommended to store anatomically shaped pads in a dry and well ventilated warehouse, at stable temperature between 5°C and 35°C and humidity level between 25% and 75%, in their original packaging. The packaging must remain intact without openings in stretch film in order to avoid dust and/or other contaminations.

## VII. OTHER

Medical device marked with CE, class I.

**Validity: 60 months**

**Issued by:** Monika Landratkiewicz, R&D Specialist

*Monika Landratkiewicz*  
  
**Technolog**

**Approved by:** Tomasz Przybylski, Deputy Director of Production and Innovation of TZMO Group

Date of issue: 14.01.2021